

NEW Micronic SampleScan Mini

A smaller, faster 2D barcode-scanning solution.



Features and Benefits:


- Scans one SBS rack in 5 - 7 seconds
- Compatible with SBS-formatted racks in 12, 24, 48, and 96 well configurations
- Built-in auto rescan feature
- Single tube read function — copies scanned barcode to the clipboard to paste into another application
- Manually edit barcodes
- Customizable output files and destination folder (.csv)
- LIMS-integration ready
- ActiveX SDK controls available to assist with robotics and LIMS integration

Software Features:


- Scans and decodes 2D barcodes
- Saves decoded data and exports output
- Excel interface
- Manual 2D barcode editing
- 384 tube compatibility option
- Auto rescan features

System Requirements:

- Windows XP computer
- 2 USB ports
- 512 MB RAM

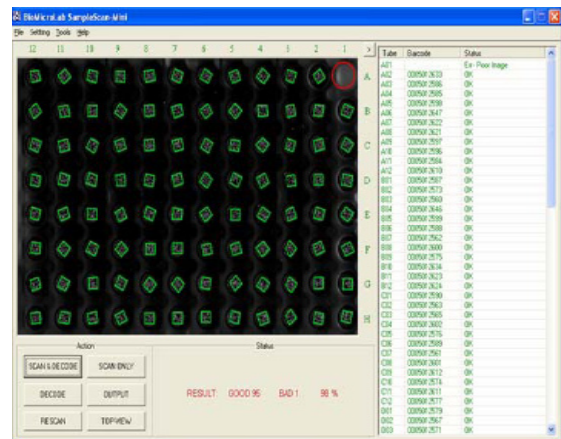

A scanned image of a rack of 96 tubes (green boxes reflect a successful code capture — red circles reflect an unsuccessful capture)

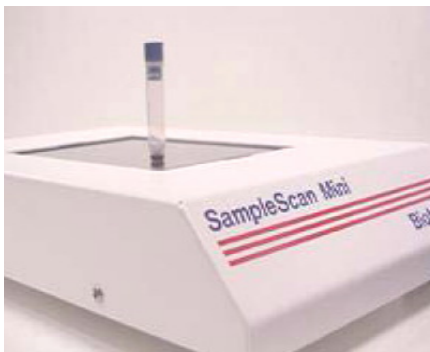


SampleScan Mini can decode racks of 96 tubes or single tubes by placing the code on the scanner surface

Specifications:

- Dimensions: 10.5" D x 6.25" W x 1.5" H
- Weight: 5 lbs.
- Electrical: 100-240 VAC, 50/60 Hz
- Resolution: 600dpi
- Warranty: One-year factory depot

SampleScan Mini is the newest addition to the Micronic line of barcode readers!

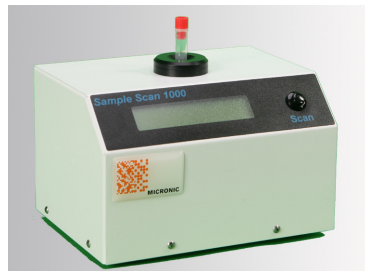

1768-5117 SampleScan Single Tube Reader



1768-5120 HandHeld Wireless Barcode Reader



1765-5104 Tracker TRxD Single Tube Reader



1768-7000 SampleScan Mini



1768-5115 SampleScan Three Rack Scanner



1768-5118 SampleScan High Speed Reader

SampleScan Mini vs. the Micronic Line of Barcode Readers

Item	Number of Racks/Tubes	384 Tube Decoding	Speed per Second	Manual Edit Mode	Capture Method	Rescan	Active X Availability	Rack Orientation Sensor	USB Ports Needed
1768-5117 SampleScan 1000	1 tube	1	Yes	Yes	Camera	No	No	N/A	3
1768-7000 SampleScan Mini	1 rack	5-7	Yes	Yes	Scanner	Yes	Yes	No	2
1768-5115 SampleScan 96	1-3 racks	20-30	Yes	Yes	Scanner	Yes	Yes	No	2
1768-5118 SampleScan HS	1 rack	3-5	Yes	Yes	Camera	Yes	Yes	Yes	3
1768-5120 Wireless HandHeld Barcode Reader	1 tube	1	No	No	Scanner	No	No	No	1
1765-5104 Tracker TRxD	1 tube	1	No	No	Camera	No	No	No	1



Please see www.usascientific.com for more information on Micronic tubes, capping options, scanners, and tube sorters.